

*ENDO*FIX^{exo}

MANUAL SUPPORT ARM





THE EASY TO USE
SUPPORT ARM FOR
GUIDING SCOPES



*ENDOFIX*exo

Steady field of vision

Precise scope guidance

Quick and easy set up

HYBRID LAPAROSCOPY

**100% German
Engineering**

ENT Surgery

STABLE VISION

THORACIC SURGERY

Throat

MANUAL CONTROL

20 years of experience

SCULLBASE SURGERY

TREMOR-FREE CAMERA GUIDANCE SYSTEM

Open Surgery

NOSE & EAR

FIELDS OF APPLICATION

ENDOFIX^{exo}




ENT Surgery

Tympanoplasty
Endoscopic Laryngeal Microdissection
Functional Endoscopic Sinus Surgery (FESS)
Endoscopic Thyroplasty
Endoscopic Removal of Nasal Polyps
Stapedectomy
Endoscopic Turbinate Reduction
Septoplasty
Endoscopic Vocal Cord Injection



Skullbase Surgery

Endoscopic Pituitary Tumor Removal /
Transsphenoidal Surgery
Endoscopic Removal of Sinonasal Tumor
Endoscopic Resection of Anterior Skull Base Tumor

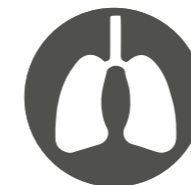

Thoracic Surgery

Wedge Resection
Esophagectomy
Lobectomy



Open Surgery

The ENDOFIX^{exo} can be used wherever the use
of a scope is indicated.

BENEFITS FOR THE PATIENT



Positive effect on surgical result

Using the ENDOFIX^{exo} enables the surgeon to work in an ergonomic and comfortable way. The surgeon can operate with both hands in a relaxed manner, which has a positive effect on the surgical result.^{1,3}

Less postoperative risk

With each 30 minutes operation time the risk of PONV (Postoperative nausea and vomiting) increases by 60 %.⁶ The patient benefits from a short operating time.

Reduced waiting times

The ENDOFIX^{exo} assistance is always available. As our system can be used for procedures without an assistant, waiting times due to staff shortages can be significantly reduced.^{1,2}

BENEFITS FOR THE CLINIC



Better scheduling efficiency

Thanks to the ENDOFIX^{exo} the surgeon guides the scope him/ herself and simple procedures can be performed without further assistance.^{1,2}

Since the ENDOFIX^{exo} is always available, it is the ideal partner in case of emergencies or staff shortages.¹

Maintenance-free and compatible

The ENDOFIX^{exo} is maintenance-free. It is also compatible with most commercially available operating tables and endoscopes.

More efficiency in OR

As efficient working is possible, valuable operating time is saved.²

Easy transport and storage

The ENDOFIX^{exo} TROLLEY allows easy transport between OR and space-saving storage when not in use.

Enhancing your clinic's image

A modern and positive image is one of the most important competitive advantages. Show your progressiveness by using robotic surgery.

BENEFITS FOR THE SURGEON



90 % of the surveyed surgeons report a significant relief in the OR thanks to robotic assistance.³

Stable positioning, steady image

Thanks to its design, the ENDOFIX^{exo} is stable and vibration-free. It does not shake or tilt the horizon, the image is steady, enabling fatigue-free work.^{1,3}

Better vision

Guiding the scope using the CONTROL is very sensitive and precise. This minimizes soiling of the front lens.

Ergonomic and relaxed working

Tensions in the neck and back, eye fatigue and the resulting pain and nausea can be significantly reduced.⁵

The ENDOFIX^{exo} supports the endoscope from above, leaving ample room for your instruments. There is no need to alter your workflow.^{1,3}

Straightforward handling

Working with the ENDOFIX^{exo} is very easy and intuitive. It is quickly attached on the operating table and covered with a STERILE COVER while the patient is being prepared for surgery.^{1,2}

Improved process coordination

As the surgeon positions the scope himself, there are no misunderstandings. The assistant can still take over the camera guidance for learning purposes.^{1,2}

Bimanual techniques welcome

The ENDOFIX^{exo} brings all relevant devices within reach of the surgeon's hands, allowing the surgeon to use one- and/ or two-handed techniques, which has a positive effect on the surgical result.¹

1. Hintschich, C.A., Fischer, R., Seebauer, C. et al. A third hand to the surgeon: the use of an endoscope holding arm in endonasal sinus surgery and well beyond. Eur Arch Otorhinolaryngol 279, 1891–1898 (2022). <https://doi.org/10.1007/s00405-021-06935-x>

2. Giotakis AI, Giotakis EI, Kyrodimos E. The Value of the Endoscope-Holding Arm in Transoral Pharyngeal Surgery. Journal of Clinical Medicine. 2024; 13(2):507. <https://doi.org/10.3390/jcm13020507>

3. Park, J.-O., Kim, M.R., Park, Y.J., et al. Transoral endoscopic thyroid surgery using robotic scope holder: Our initial experiences. J Minim Access Surg. 2020 Jul-Sep;16(3):235-238. doi: 10.4103/jmas.JMAS_12_19. PMID: 31031326; PMCID: PMC7440021

4. Lin YY, Hsieh MJ, Wu CY, Yang LY, Pan YB, Wu CF, Gonzalez-Rivas D, Chao YK. Comparison of active versus passive robotic-endoscope-holder-assisted unisurgeon uniportal thoracoscopic surgery in terms of surgical efficacy and patient safety. J Thorac Dis 2023;15(7):3800-3810. doi: 10.21037/jtd-23-19

5. Dixon F, Vitish-Sharma P, Khanna A, Keeler BD, Work-related musculoskeletal pain and discomfort in laparoscopic surgeons: an international multispecialty survey. Ann R Coll Surg Engl 2023; 105: 734–738. <https://doi.org/10.1308/rcsann.2023.0024>

6. Sinclair DR, Chung F, Mezei G: Can Postoperative Nausea and Vomiting Be Predicted? Anesthesiology. 1999;91(1):109–18. Doi:10.1097/0000542-199907000-00018



ABOUT OUR SYSTEM



ENDOFIX^{exo}

Manual, easy to use support arm for holding and guiding scopes

In modern ENT and Skullbase Surgery, microscopic techniques are increasingly replaced by endoscopic surgery. When working endoscopically, the surgeon's field of vision increases significantly* and his/ her position becomes more ergonomic.

With the development of the ENDOFIX ^{exo}, an intuitively movable manual support arm is now available. It ensures the stable and sensitive positioning of the scope and enables the surgeon to use one and/ or two hands techniques, which has a positive effect on the surgical result.

Since the ENDOFIX ^{exo} can be used with both, endoscope and exoscope, it is also the perfect assistant in Open Surgery.

ABOUT OUR SYSTEM



Optimum usability

In minimally invasive procedures, particularly in ENT and Skullbase Surgery, the surgeon prefers to work alone and with his/ her two hands. We have developed the ENDOFIX exo to perfectly meet these requirements.


The ENDOFIX exo is equipped with a counter-balance to compensate the weight of the scope and the camera head.

Complete weight compensation is achieved in a middle cantilever position. The scope is held seemingly weightless.


The ENDOFIX exo is supporting the scope from above leaving the surgeon sufficient space for instruments and ergonomic two-handed working. An ideal condition for comfortable surgical work.

ABOUT OUR SYSTEM

Mount setting for scopes with cylindrical eyepiece Ø7-20 mm

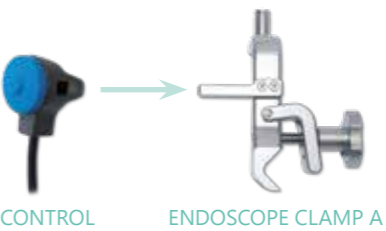

Mount setting for scopes with Ø10 mm or Storz-type cone

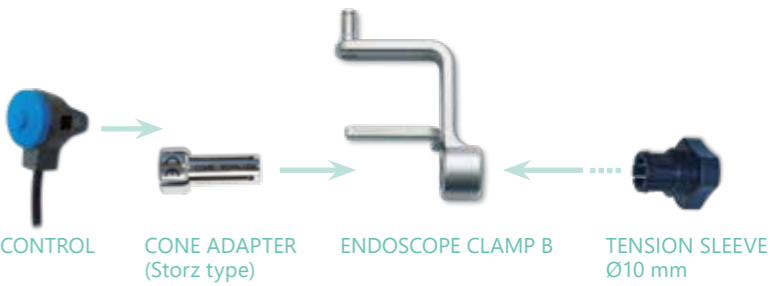

How to control the ENDOFIX exo ?

The ENDOFIX exo is ready to operate once the green LED on the DEACTIVATION PUSHBUTTON is lit. Before surgery, the arm is moved and pulled into the required position by pressing the DEACTIVATION PUSHBUTTON. When the button is released, the arm is locked and remains in the set position. During surgery the ENDOFIX exo is operated with the CONTROL.

The CONTROL is designed for intuitive and ergonomic operation of the ENDOFIX exo. Thanks to the CONTROL the scope can be moved very precisely. The CONTROL is connected to the ENDOSCOPE CLAMPS within reach of the surgeon's fingers. Both, CONTROL and ENDOSCOPE CLAMPS are designed to be reprocessed.

ABOUT THE ENDOFIX^{exo}



A finely tuned system of components allows the adaptation to most commercially available scopes. The CONTROL and the scopes can be attached to the ENDOFIX^{exo} with the ENDOSCOPE CLAMPS.


The components of the ENDOFIX^{exo} like CONTROL, ENDOSCOPE CLAMP A, ENDOSCOPE CLAMP B, TENSION SLEEVE and CONE ADAPTER are multi-use and can be reprocessed in our STERILISATION TRAY.

Before every intervention the ENDOFIX^{exo} is covered by a single-use STERILE COVER.

After the surgical procedure, the ENDOFIX^{exo} is hung onto the TROLLEY which is designed for easy and safe transport and storage of our systems.



Quick-coupling device for easy mounting to any surgical table with standard rail

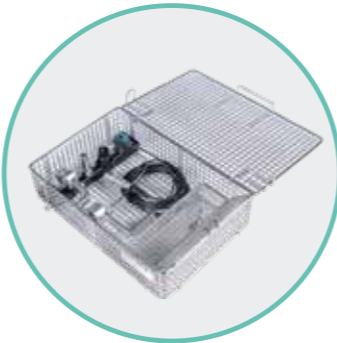


STERILISATION TRAY for reprocessing our autoclavable components

SYSTEM COMPONENTS & ACCESSORIES

	ENDOFIX^{exo} Support Arm / TROLLEY Thanks to its quick-coupling device for easy mounting to any surgical table with standard rail. The TROLLEY is designed for safe transport and storage of our robotic arms.
	CONTROL * For controlling and positioning the ENDOFIX ^{exo}
	ENDOSCOPE CLAMP A * / ENDOSCOPE CLAMP B * For attaching the scope and CONTROL to the ENDOFIX ^{exo}
	TENSION SLEEVE D10 * For Ø10 mm diameter scopes
	CONE ADAPTER * For scopes with Storz type cone
	STERILISATION TRAY * Developed for simultaneous autoclaving of all reprocessable components
	STERILE COVER Before every intervention the ENDOFIX ^{exo} is covered by a single-use STERILE COVER.

* autoclavable part


ORDERING THE ENDOFIX^{exo} SYSTEM*

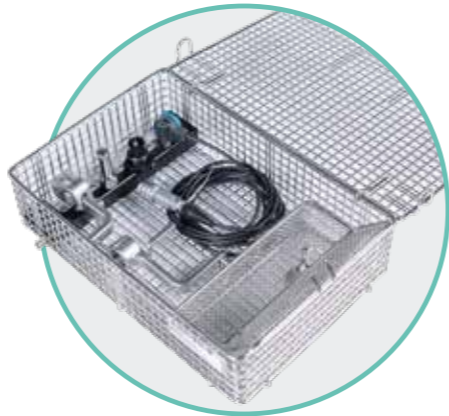

Article number	Designation	Content / Description
171912	ENDOFIX EXO SET	<ul style="list-style-type: none">• 1 ENDOFIX^{exo}• 1 TROLLEY
171909	BUNDLE EXO CONTROL	<ul style="list-style-type: none">• 1 CONTROL• 3 STERILE COVER PU/ 50 ea.
182233	ENDOSCOPE CLAMP A	For attaching scopes with cylindrical eyepiece Ø7-20 mm
151686	ENDOSCOPE CLAMP B	For attaching scopes with Ø10 mm diameter or Storz-type cone
141339	TENSION SLEEVE D10	For connecting scopes with Ø10 mm diameter and for CONE ADAPTER; with ENDOSCOPE CLAMP B
171927	CONE ADAPTER	For attaching scopes with Storz-type cone to TENSION SLEEVE D10
110084	STERILE COVER PU/ 50 ea.	Pack of 50 individually packed STERILE COVER for ENDOFIX ^{exo}
141527	STERILISATION TRAY	Special STERILISATION TRAY for ENDOFIX ^{exo} sterilisable parts including a SMALL PARTS TRAY
171906	CONTROL	For controlling and positioning the ENDOFIX ^{exo}


TECHNICAL SPECIFICATION

All dimensions ± 5 mm

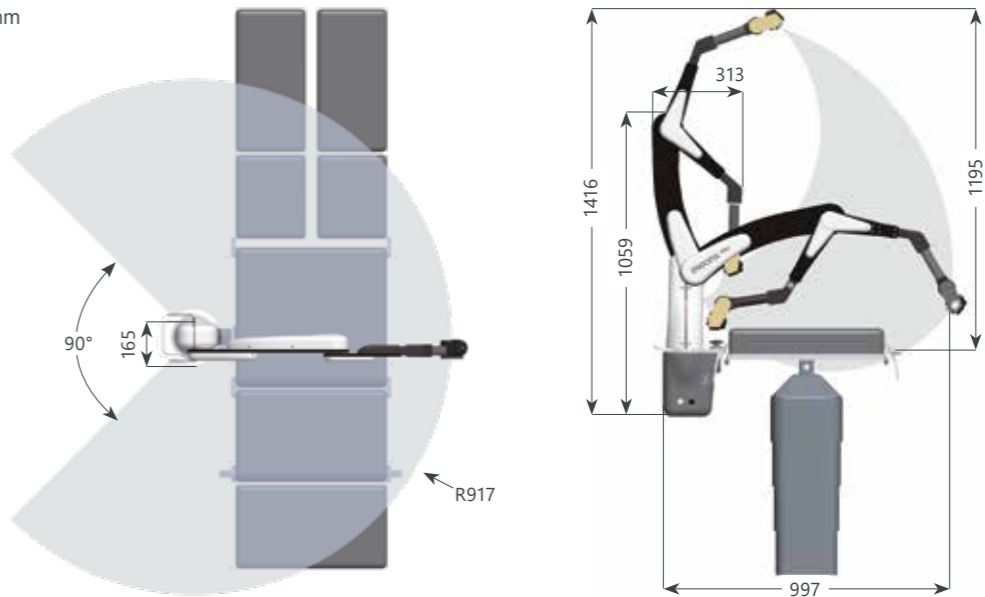

Feature	ENDOFIX ^{exo}
Approval	CE (class I), FDA 510 (k)
Supply voltage	100 - 240 Volt, 50 - 60 Hz
Maximum power consumption	60 VA
Size	313 x 165 x 1059 mm
Weight	12.0 kg
Safe working load	1.0 kg
Lockable axis	6
CONTROL & ENDOSCOPE CLAMPS	Multiuse
Endoscope size	7 to 20 mm cylindrical eyepiece, 10 mm or with Storz-type cone
OR table adaption	Quick fastener, suitable for OR tables with European and US size rails, UK rails with additional adapters

* This list is only an excerpt and not exhaustive.

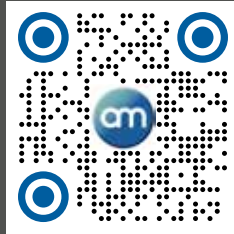

Developed & Manufactured by

AKTORMed GmbH
Neugablonzer Strasse 13
93073 Neutraubling, Germany

Email: sales@aktormed.com
www.aktormed.com

