

POSITION GUIDE

FOR ROBOTIC CAMERA CONTROL

SOLOASSIST^{II}



INTRODUCTION

Laparoscopic techniques have become the standard for many procedures. Perfect imaging is the basis for safe and effective surgical action. Therefore the stable holding and flexible positioning of the endoscope is of great importance.

The SOLOASSIST supports this in a unique way. In order to obtain the best possible result, it is necessary to position the SOLOASSIST on the operating table in such a way that:

- The situs can be completely viewed.
- Neither the surgeon, nor assistants are obstructed.
- The direct view to the monitor is not hidden.
- No collisions occur between the SOLOASSIST and the surgeon's instruments.

Thanks to the high flexibility of the SOLOASSIST, in most cases it is possible to guarantee the above mentioned points.

This document helps you to find the ideal position of the SOLOASSIST for different procedures. The shown examples are only suggestions. In practice, every user has a slightly different setup, e.g. regarding the trocar positions.

Before the first operation, the positioning of the SOLOASSIST should be discussed with the surgeon. The aim is not for the surgeon to adapt to the device, but for the SOLOASSIST to be positioned in such a way that the surgeon can work with his usual procedure.

ATTENTION! The trocar positions shown in this guide are for illustration purposes only and do not represent a recommendation. The positioning of the trocars is the sole responsibility of the surgeon.



AGENDA

SYSTEM OVERVIEW	4
BRIEF SETUP INSTRUCTION SOLOASSIST II	5
SOLOASSIST II POSITIONING	6
LAPAROSCOPIC APPENDECTOMY	8
LAPAROSCOPIC INGUINAL HERNIA REPAIR	9
LAPAROSCOPIC INGUINAL HERNIA REPAIR (BILATERAL)	10
LAPAROSCOPIC CHOLECYSTECTOMY	11
LAPAROSCOPIC CHOLECYSTECTOMY (ALTERNATIVE)	12
LAPAROSCOPIC CHOLECYSTECTOMY (SINGLE INCISION)	13
LAPAROSCOPIC DISTAL GASTRECTOMY	14
LAPAROSCOPIC TOTAL GASTRECTOMY	15
LAPAROSCOPIC RIGHT COLECTOMY	16
LAPAROSCOPIC LEFT COLECTOMY	17
LAPAROSCOPIC RECTAL RESECTION	18
LAPAROSCOPIC DISTAL PANCREATECTOMY	19
LAPAROSCOPIC FUNDOPLICATION	20
LAPAROSCOPIC SIGMA RESECTION	21
LAPAROSCOPIC GASTRIC SLEEVE	22
LAPAROSCOPIC PROSTATECTOMY	23
LAPAROSCOPIC NEPHRECTOMY RIGHT	24
LAPAROSCOPIC NEPHRECTOMY RIGHT (ALTERNATIVE)	25
LAPAROSCOPIC HYSTERECTOMY	26
LAPAROSCOPIC OVARIAN CYSTECTOMY	27

SYSTEM OVERVIEW



BRIEF SETUP INSTRUCTION

Fasten the SOLOASSIST II to the OP table.




Connect the cable of the control and the power supply to the jacks.



Cover the robotic arm with a STERILE COVER.



Save the TROCAR POINT.
Press the UNLOCK BUTTON, hold it while you lead the TRACER PIN to the point of entry of the trocar. Then press the button "TP" to save the TROCAR POINT.

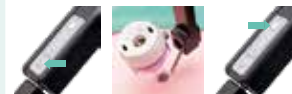

Mount the ENDOSCOPE CLAMP and the TENSION SLEEVE to the endoscope.



Attach the endoscope to the SOLOASSISTII.



Attach the JOYSTICK to a surgical instrument.




PLEASE NOTE: THIS BRIEF INSTRUCTION DOES NOT REPLACE THE CAREFUL READING OF THE OPERATION MANUAL.

CAUTION: ONLY INTRODUCED STAFF IS AUTHORIZED TO HANDLE THE SOLOASSISTII.

TRY TO INDICATE THE FOLLOWING

- Position of patient
- Position of assistant and surgeon
- Position of patient arms (adjacent, extended)
- Position of first and second assistant
- Camera port ●
- Operation port position ●
- Support port position ●
- Target organ ★
- Surgeons instrument →
- Endoscope →
- Assistants instrument →

CW: CLOCKWISE
CCW: COUNTER CLOCKWISE
SAIL: SOLOASSIST II

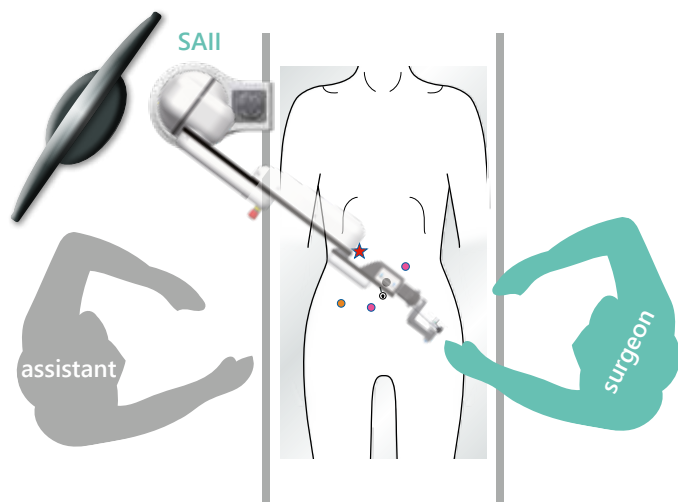


HOW TO GET THE RIGHT POSITION

The following basic rules makes it easy to find a suitable position for the SOLOASSIST. This position may be optimized according to the experience gained during the first interventions.

- Place the SOLOASSIST opposite to the surgeon in a position that a line is formed from the surgeon via the trocar to the device.
- Choose joint position to get a right angle between scope and universal joint.


How to change the joint position

1. Pull the joint head forward.
2. Turn the pivot joint into the desired position.
3. Let the joint head lock in again.

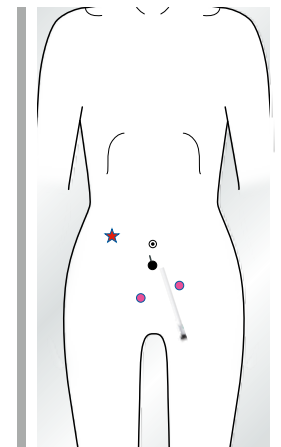


LAPAROSCOPIC APPENDECTOMY

Body position	Supine	Assistant position	-
Number of ports	3	SAll position	Lower right
Scope port	Lower left	Universal joint	0°
Surgeon position	Left	Personal note	

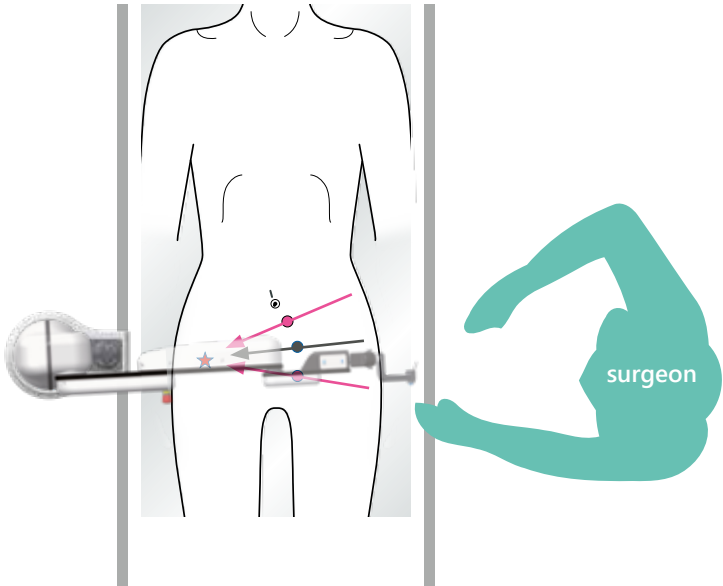


LAPAROSCOPIC INGUINAL HERNIA REPAIR

Body position	Supine	Assistant position	-
Number of ports	3	SAll position	Lower opposite side
Scope port	Umbilicus	Universal joint	45° CW right side hernia 45° CCW left side hernia
Surgeon position	Opposite to affected side	Personal note	

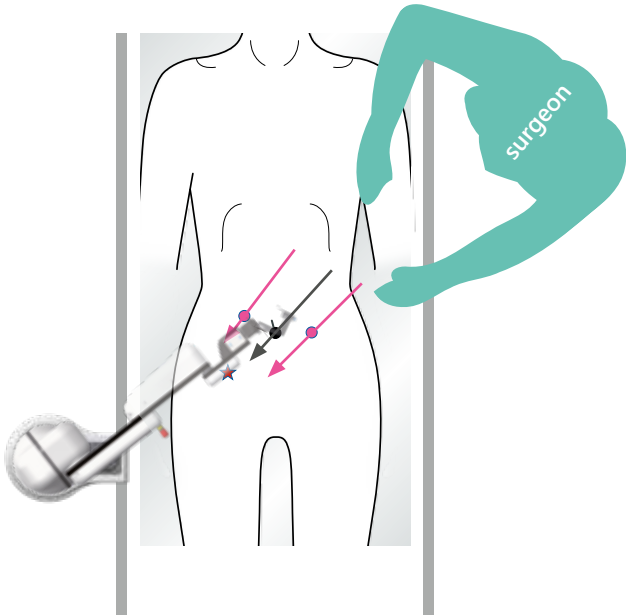


LAPAROSCOPIC INGUINAL HERNIA REPAIR (BILATERAL)

Body position	Supine	Assistant position	-
Number of ports	3	SAll position	Lower opposite to first side
Scope port	Umbilicus	Universal joint	45° CW
Surgeon position	Opposite to first side	Personal note	

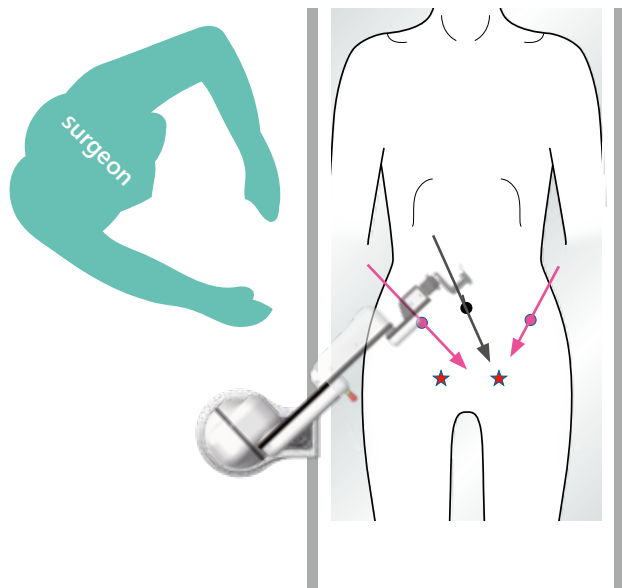


LAPAROSCOPIC CHOLECYSTECTOMY

Body position	Supine	Assistant position	Right
Number of ports	4	SAll position	Upper right
Scope port	Umbilicus	Universal joint	45° CW
Surgeon position	Left	Personal note	

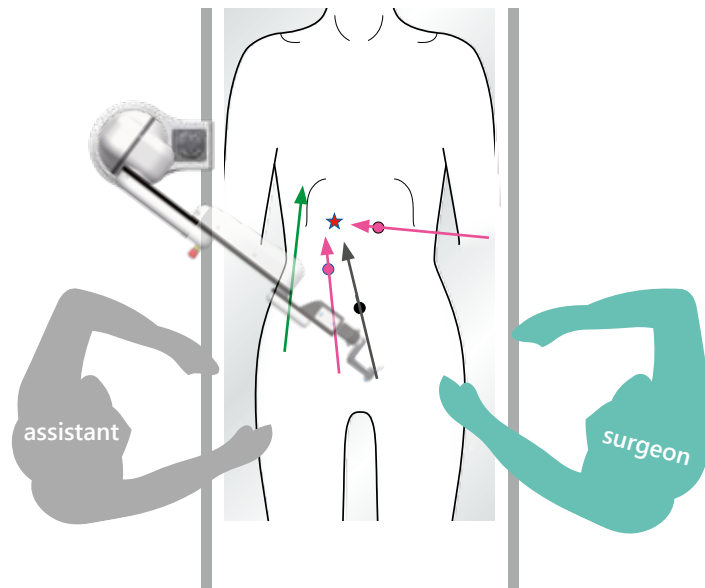


LAPAROSCOPIC CHOLECYSTECTOMY (ALTERNATIVE)

Body position	Supine	Assistant position	Right
Number of ports	4	SAll position	Lower right
Scope port	Umbilicus	Universal joint	90° CW
Surgeon position	Left	Personal note	

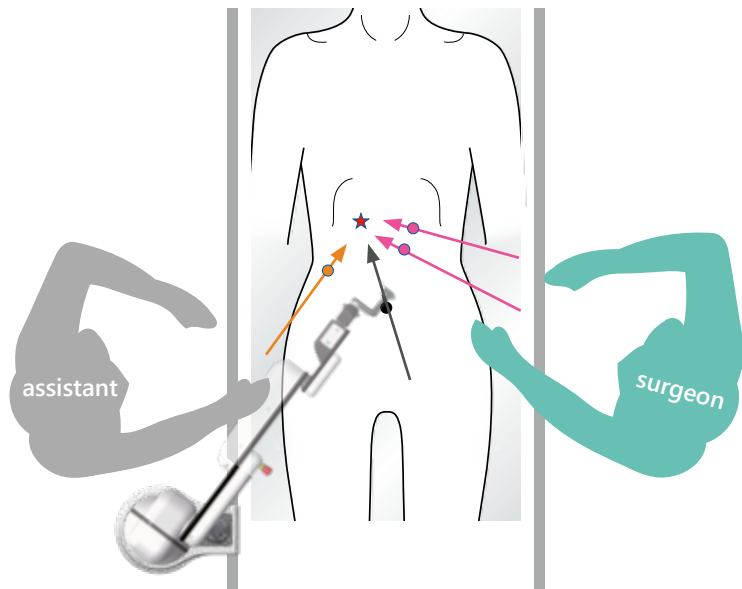


LAPAROSCOPIC CHOLECYSTECTOMY (SINGLE INCISION)

Body position	Supine	Assistant position	-
Number of ports	1	SAll position	Upper left
Scope port	Umbilicus	Universal joint	135° CW
Surgeon position	Left	Personal note	

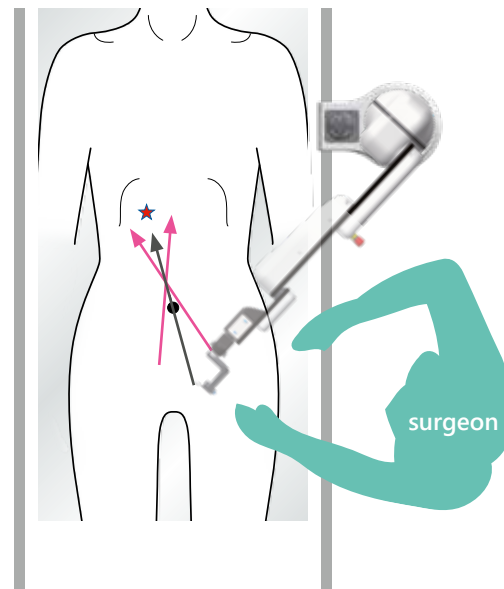


LAPAROSCOPIC DISTAL GASTRECTOMY

Body position	Lithotomy	Assistant position	Right
Number of ports	4	SAll position	Upper left
Scope port	Umbilicus	Universal joint	45° CCW
Surgeon position	Between the legs	Personal note	

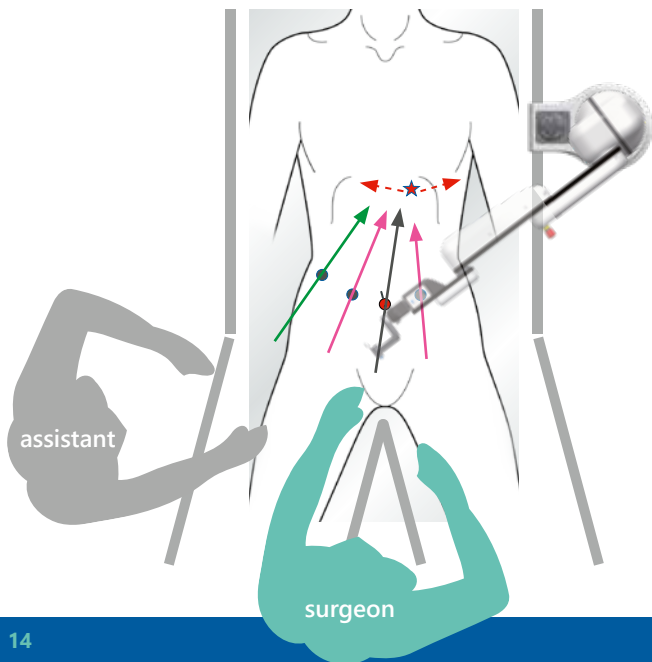


LAPAROSCOPIC TOTAL GASTRECTOMY

Body position	Lithotomy	Assistant position	Right
Number of ports	4	SAll position	Upper left
Scope port	Umbilicus	Universal joint	45° CCW
Surgeon position	Between the legs	Personal note	

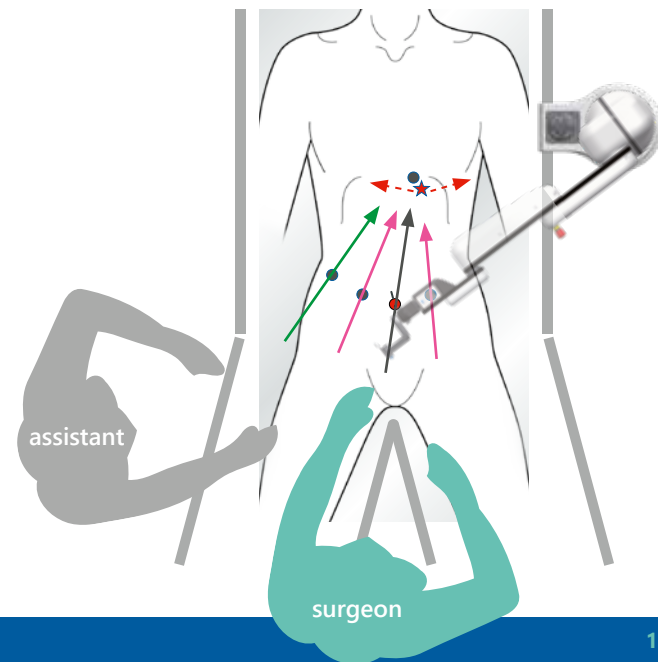


LAPAROSCOPIC RIGHT COLECTOMY

Body position	Supine	Assistant position	Right
Number of ports	4	SAll position	Upper right
Scope port	Umbilicus	Universal joint	0°
Surgeon position	Left	Personal note	

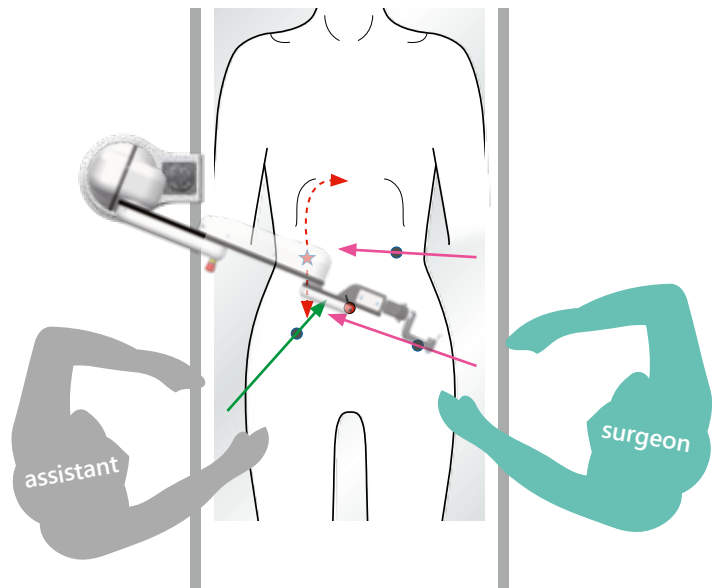


LAPAROSCOPIC LEFT COLECTOMY

Body position	Supine	Assistant position	Left
Number of ports	4	SAll position	Lower left
Scope port	Umbilicus	Universal joint	0°
Surgeon position	Right	Personal note	

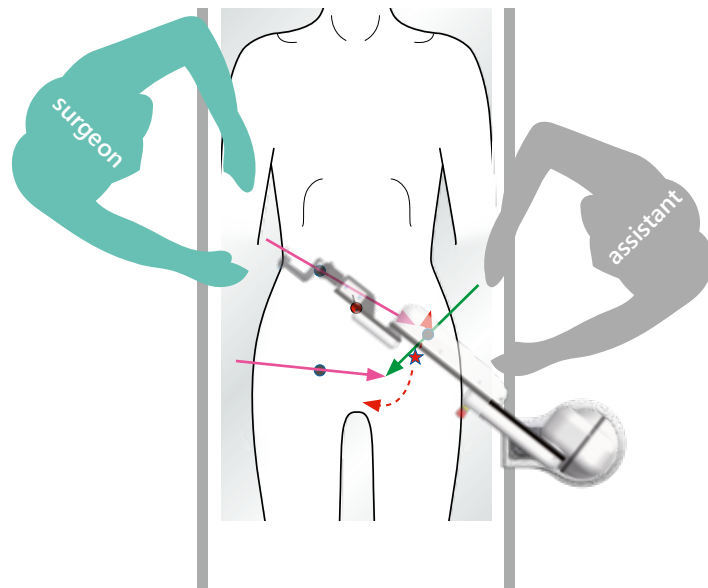


LAPAROSCOPIC RECTAL RESECTION

Body position	Supine	Assistant position	Left
Number of ports	5	SAll position	Upper left
Scope port	Umbilicus	Universal joint	45° CW
Surgeon position	Right	Personal note	

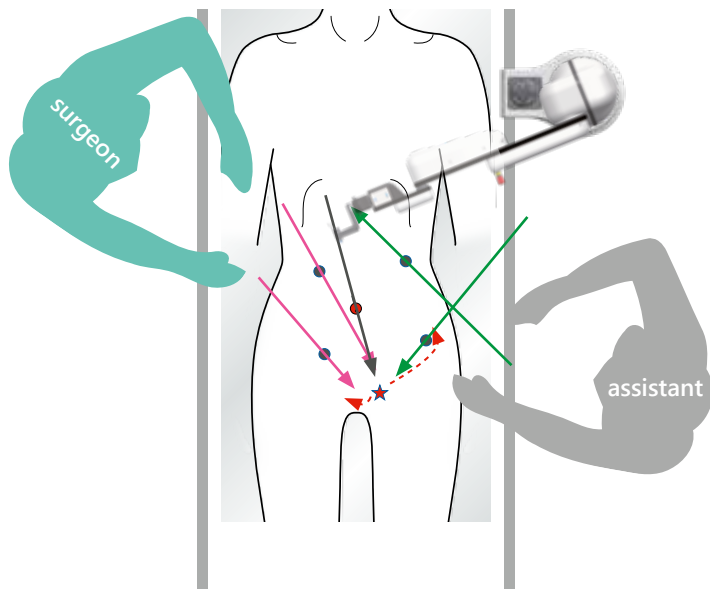


LAPAROSCOPIC DISTAL PANCREATECTOMY

Body position	Lithotomy	Assistant position	Right
Number of ports	4	SAll position	Upper left
Scope port	Umbilicus	Universal joint	45° CCW
Surgeon position	Between the legs	Personal note	

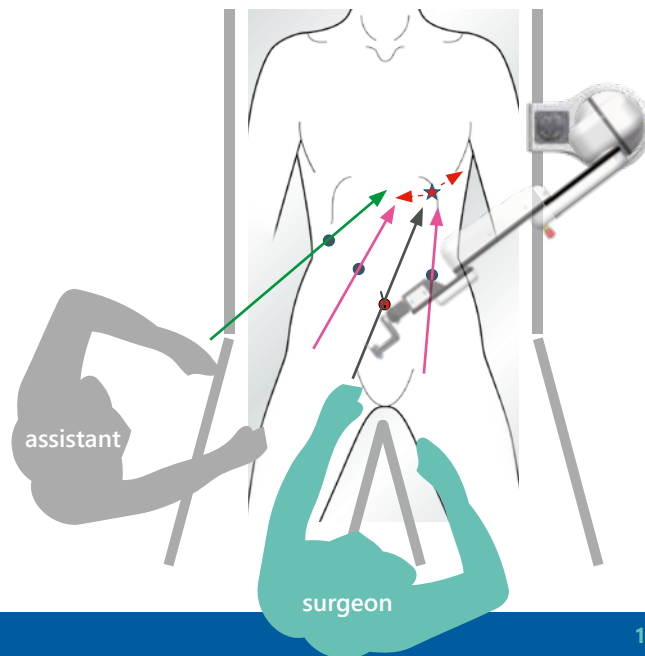


LAPAROSCOPIC FUNDOPLICATION

Body position	Lithotomy	Assistant position	-
Number of ports	3 (4)	SAll position	Left
Scope port	Umbilicus	Universal joint	45° CW
Surgeon position	Between the legs	Personal note	

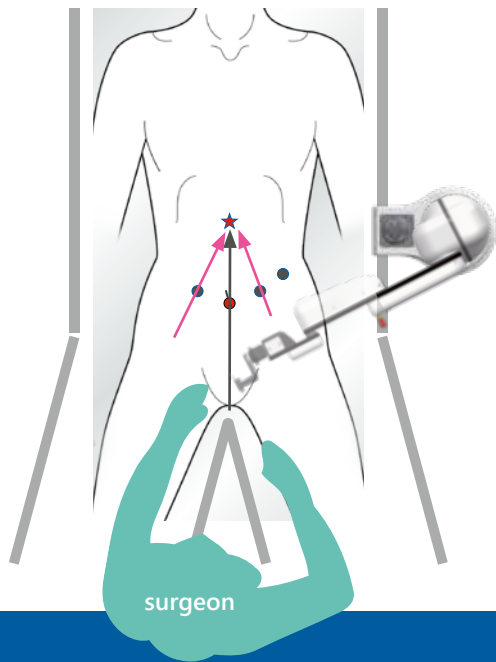


LAPAROSCOPIC SIGMA RESECTION

Body position	Lithotomy	Assistant position	-
Number of ports	3 (4)	SAll position	Left
Scope port	Umbilicus	Universal joint	0°
Surgeon position	Right	Personal note	

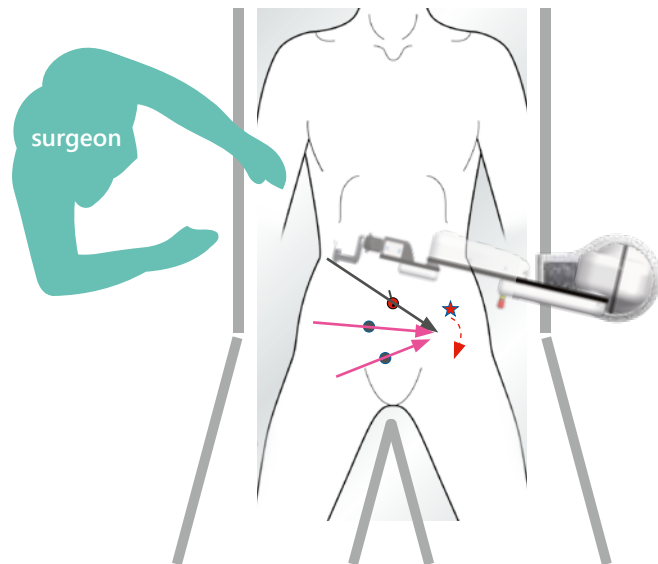


LAPAROSCOPIC GASTRIC SLEEVE

Body position	Lithotomy	Assistant position	Right
Number of ports	4	SAll position	Upper left
Scope port	Umbilicus	Universal joint	45° CCW
Surgeon position	Between the legs	Personal note	

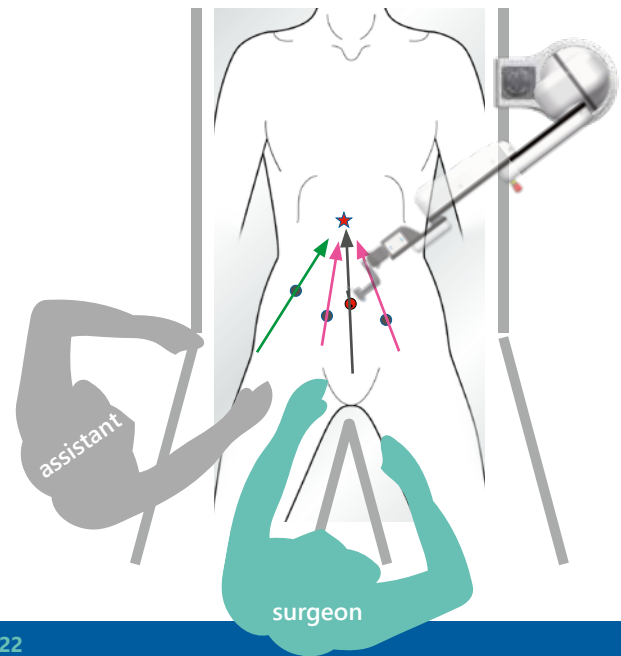


LAPAROSCOPIC PROSTATECTOMY

Body position	Supine	Assistant position	Right
Number of ports	5	SAll position	Lower right
Scope port	Umbilicus	Universal joint	45° CW
Surgeon position	Left	Personal note	

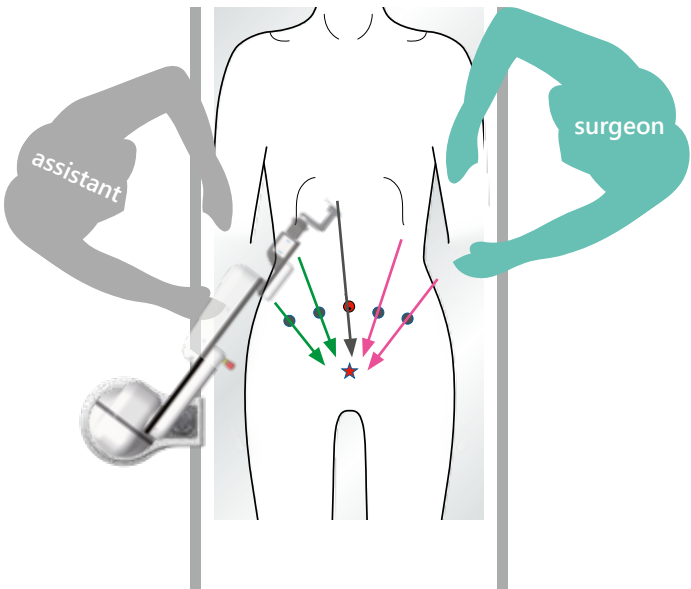


LAPAROSCOPIC NEPHRECTOMY RIGHT

Body position	45° Side	Assistant position	-
Number of ports	3	SAll position	Right
Scope port	Umbilicus	Universal joint	0°
Surgeon position	Left	Personal note	

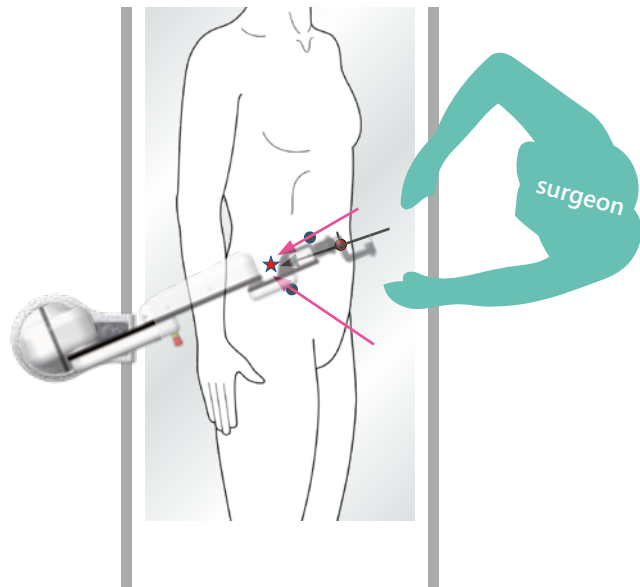


LAPAROSCOPIC NEPHRECTOMY RIGHT (ALTERNATIVE)

Body position	45° Side	Assistant position	-
Number of ports	3	SAll position	Upper left
Scope port	Umbilicus	Universal joint	180° CW
Surgeon position	Left	Personal note	

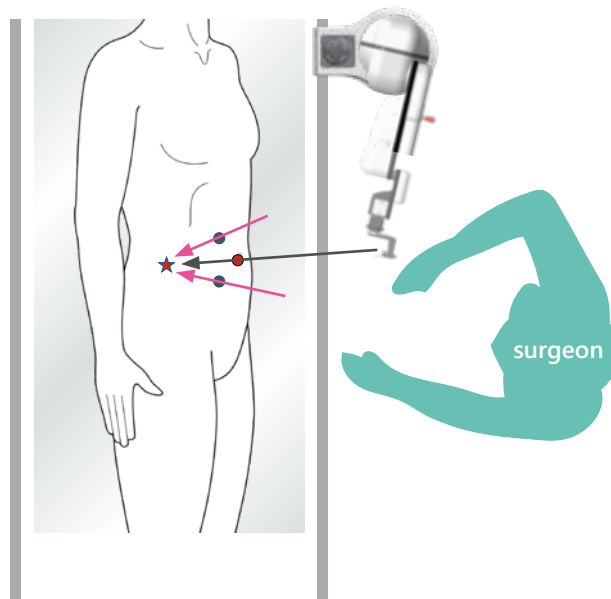


LAPAROSCOPIC HYSTERECTOMY

Body position	Lithotomy	Assistant position	Right
Number of ports	4	SAll position	Upper right
Scope port	Umbilicus	Universal joint	45° CCW
Surgeon position	Left	Personal note	

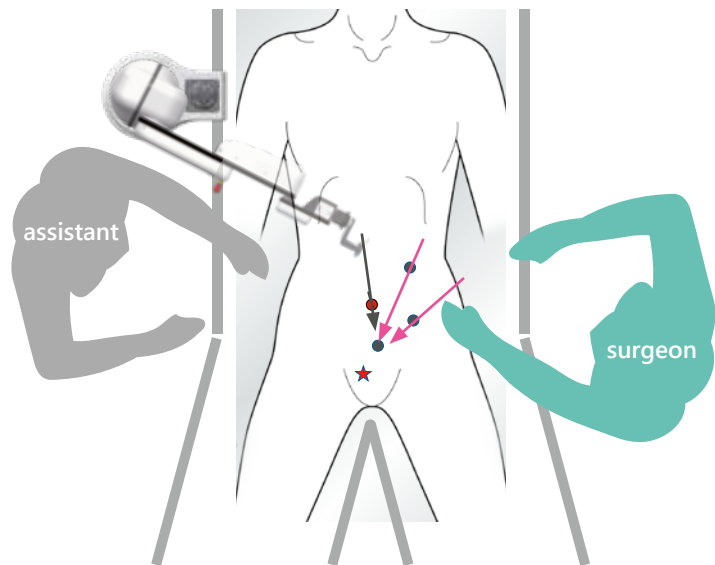


LAPAROSCOPIC OVARIAN CYSTECTOMY

Body position	Supine	Assistant position	-
Number of ports	3	SAll position	Lower right
Scope port	Umbilicus	Universal joint	45° CCW
Surgeon position	Upper left	Personal note	

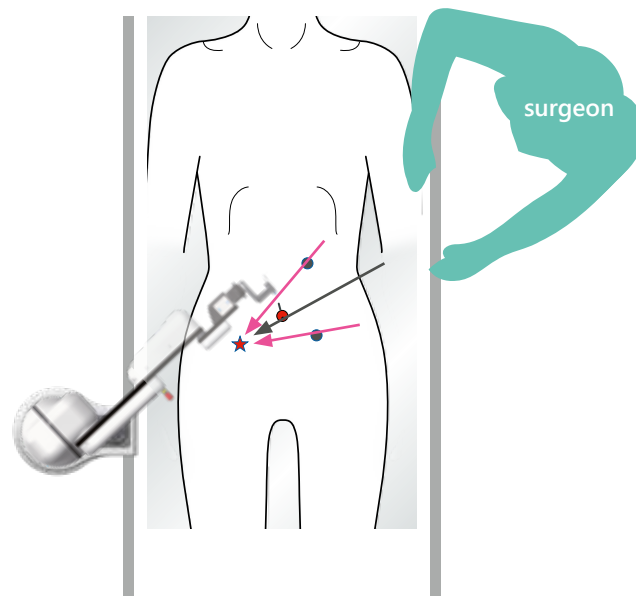


YOUR PROCEDURE

Body position		Assistant position	
Number of ports		SAll position	
Scope port		Universal joint	
Surgeon position		Personal note	

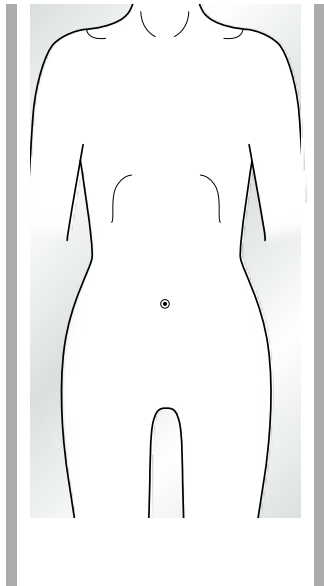


YOUR PROCEDURE

Body position		Assistant position	
Number of ports		SAll position	
Scope port		Universal joint	
Surgeon position		Personal note	

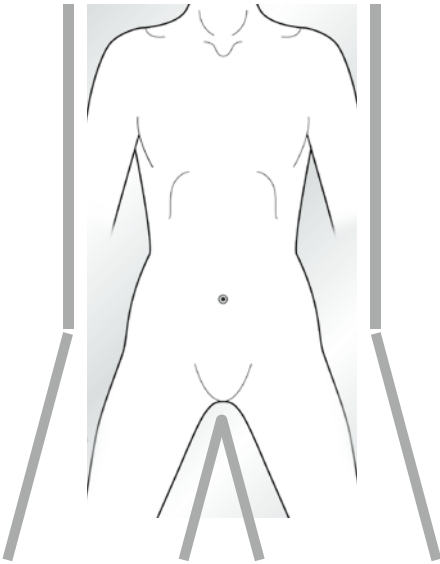


YOUR PROCEDURE

Body position		Assistant position	
Number of ports		SAll position	
Scope port		Universal joint	
Surgeon position		Personal note	

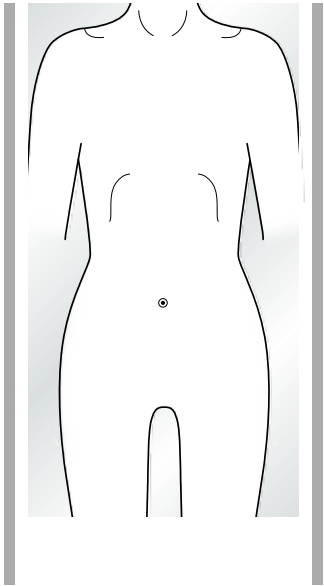

YOUR PROCEDURE

Body position		Assistant position	
Number of ports		SAll position	
Scope port		Universal joint	
Surgeon position		Personal note	

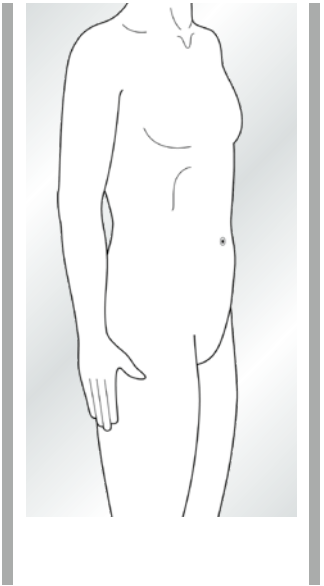

NYX, managed document! Number: 0103-183618 Rev: 04, freigegeben / released: Robert Geiger 22.09.2025
Printed versions of this document are for information only!



Developed & Manufactured by

AKTORMed GmbH
Neugablonzer Strasse 13
93073 Neutraubling, Germany

Email: sales@aktormed.com

www.aktormed.com



This brochure is not for distribution in the US.

0103-183618-04